

Product Description

Product Name: HCT116-hNIS-Neo/iRFP-Puro
Catalog Number: CL095
Lot Number: CL-IM93
Species: Human (*Homo sapiens*)
Tissue: Colon
Cell type: Colorectal carcinoma
Parental cells: HCT116 (ATCC® CCL-247™)*
Morphology: Epithelial
Growth mode: Adherent
Reporter genes: Human sodium iodide symporter (hNIS)
Near-infrared fluorescent protein (iRFP)
Selection genes: Neomycin (Neo)
Puromycin (Puro)

This is a polyclonal population derived from the colorectal carcinoma HCT116 cell line (ATCC® CCL-247™). Parental HCT116 cells were transduced with 1) LV-hNIS-IRES-Neo (Imanis #LV013) encoding the human sodium iodide symporter (hNIS) cDNA under the spleen focus-forming virus (SFFV) promoter linked to the neomycin resistance gene (Neo) via an IRES and 2) LV-iRFP-P2A-Puro (Imanis #LV032) encoding the near-infrared fluorescent protein (iRFP; ex/em = 690/713) cDNA under the SFFV promoter linked to the puromycin resistance gene (Puro) via a P2A cleavage peptide. High hNIS and iRFP expressing cells were selected using G418 and puromycin. The lentiviral vectors are self-inactivating (SIN) vectors in which the viral enhancer and promoter have been deleted. Transcription inactivation of the LTR in the SIN provirus increases biosafety by preventing mobilization by replication competent viruses and enables regulated expression of the genes from the internal promoters without *cis*-acting effects of the LTR¹.

* The ATCC trademark and any and all ATCC catalog numbers are trademarks of the American Type Culture Collection

Mycoplasma Testing

The HCT116-hNIS-Neo/iRFP-Puro cell line has been tested for mycoplasma contamination and is certified mycoplasma free.

Cell Line Authentication

In light of studies suggesting that 18-36% of cell lines utilized in biomedical research are contaminated or completely misidentified,^{2,3} several funding organizations, including NIH, as well as major publishers, including those affiliated with the American Association for Cancer Research (AACR), require cell lines used in research to be authenticated prior to publication^{4,5}. The parental HCT116 cell line used to generate HCT116-hNIS-Neo/iRFP-Puro was authenticated and certified free of interspecies cross contamination by STR profiling with 9 STR loci.

Recommended Uses

HCT116-hNIS-Neo/iRFP-Puro cells are suitable for *in vitro* and *in vivo* experimentation. HCT116 cells form primary tumors and distant metastases (lung, liver, and lymph nodes) post implantation into immunosuppressed mice^{6,7}. The Fluc and iRFP transgenes facilitate high-resolution 3D PET/SPECT and noninvasive imaging of implanted cells, respectively.

References

- ¹Miyoshi et al. J Virol. 1998. 72:8150-8157.
- ²Hughes et al. BioTechniques 2007. 43: 575-586.
- ³Chatterjee et al. Science 2007. 315:928-931.
- ⁴<https://grants.nih.gov/grants/guide/notice-files/NOT-OD-08-017.html>
- ⁵<http://www.aacrjournals.org/site/InstrAuthors/fora.xhtml#celllineuse>
- ⁶Brattain et al. Cancer Res. 1981. 41:1751-1756.
- ⁷Yang et al. Anticancer Res. 1997. 17:3463-3468.

Storage Instructions

Remove cells from the dry ice packaging and immediately store in the vapor phase above liquid nitrogen (below -130°C).

Complete Growth Medium

Dulbecco's Modified Eagle's Medium (DMEM)
10% fetal bovine serum (FBS)
1% Penicillin/Streptomycin
0.5 mg/mL G418 (to maintain high hNIS expression)
1 µg/mL puromycin (to maintain high iRFP expression)

Thawing Instructions

1. Thaw cells by gently swirling in a 37°C water bath. To limit contamination, do not submerge the O-ring and cap.
2. When cells are ~70% thawed (less than 1 min), remove the vial and wipe down with 70% ethanol. Allow tube to dry completely.
3. In a biosafety cabinet, transfer the cells into a 15 mL conical tube containing 5 mL of pre-warmed complete growth medium without selection drugs. Centrifuge cells at ~250 x g for 3-5 min.
4. Remove supernatant and resuspend cells in 1 mL complete growth medium without selection drugs. Transfer cells to a T75 flask containing 10 mL pre-warmed complete growth medium without selection drugs.
5. Incubate the culture at 37°C with 5% CO₂. After 48 hours, replace the culture supernatant with complete growth medium containing 0.5 mg/mL G418 and 1 µg/mL puromycin. Cells should reach full confluency 3-4 days after thawing.

Subculturing Instructions

Volumes are given for a T75 flask; increase or decrease as needed. To maintain high hNIS and iRFP expression, it is recommended that cells be subcultured in the presence of 0.5 mg/mL G418 and 1 µg/mL puromycin. The cells should be passaged when they reach 90-100% confluency.

1. Remove culture medium from cells.
2. Carefully wash the cell monolayer with 5-10 mL of phosphate buffered saline.
3. Add 2 mL of 0.25% Trypsin-EDTA solution to the flask and incubate at 37°C until cells have dissociated (approx. 2-5 min).
4. Neutralize the trypsin by adding 8 mL complete growth medium, and mix by gently pipetting up and down.
5. Transfer desired portion of the cells to a fresh T75 flask. Add fresh complete growth medium to a total volume of 10 mL and return cells to 37°C/5% CO₂ incubator.

For maintenance a subcultivation ratio of 1:10 is recommended. At this ratio cells will be ready for passage every 3-4 days.

Freezing Medium

The cells can be amplified and used to generate additional frozen stocks. Frozen stocks should be preserved in a designated cryopreservation medium or in complete growth medium without selection drugs supplemented with 5-10% DMSO.

Additional Considerations

HCT116 cells can detach from tissue culture surfaces upon repeated washing. Coating culture plates with poly-D-Lysine prior to use can increase cell adherence to the plates during experiments requiring repeated washes or media changes.


Certificate of Analysis

Testing performed by Imanis Life Sciences


Test description	Result
Post thaw viable cell recovery	Pass QC
Sterility	No contamination detected
Mycoplasma	No contamination detected
Neomycin selection	Pass QC
Puromycin selection	Pass QC
NIS expression	Pass QC
Fluorescence expression	Pass QC

Morphology:

Low density, 200X

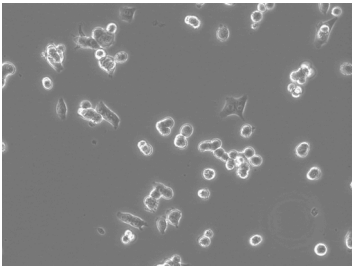


High density, 200X

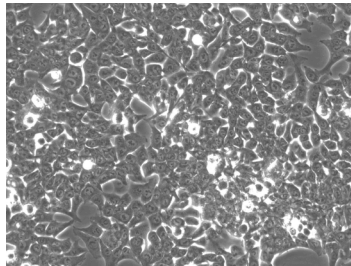


Low and high density photos taken 24 and 69 hours after thawing, respectively.

NIS Expression:

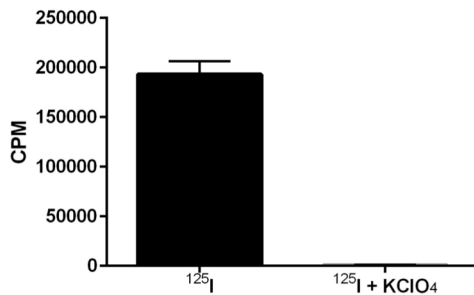

Uptake of ^{125}I by 3×10^5 cells was assayed in the presence or absence of KClO_4 , an inhibitor of NIS-mediated ^{125}I uptake.

Fluorescence Expression:

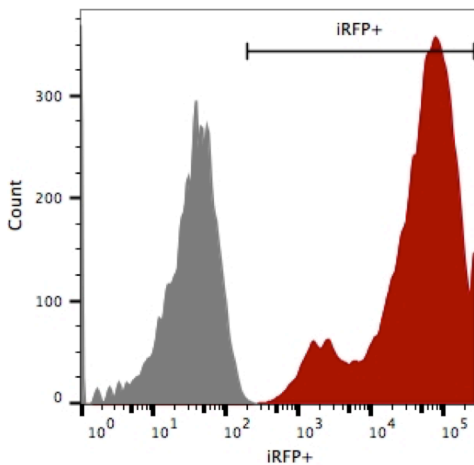

HCT116-hNIS-Neo/iRFP-Puro (red) or isotype control (HCT116; grey) cells were fixed with paraformaldehyde and analyzed by flow cytometry.

Legal Disclaimers

LIMITED PRODUCT WARRANTY

THIS WARRANTY LIMITS OUR LIABILITY TO REPLACEMENT OF THIS PRODUCT. NO OTHER WARRANTIES OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE PROVIDED BY IMANIS. IMANIS SHALL HAVE NO LIABILITY FOR ANY DIRECT, INDIRECT, CONSEQUENTIAL, OR INCIDENTAL DAMAGES ARISING OUT OF THE USE, THE RESULTS OF USE, OR THE INABILITY TO USE THIS PRODUCT.

FOR *IN VITRO* USE ONLY. THIS CERTIFICATE IS A DECLARATION OF ANALYSIS AT THE TIME OF MANUFACTURE.

PURCHASER NOTIFICATION

LIMITED LICENSE NOTICE – RESEARCH USE ONLY

IMANIS LIFE SCIENCES HAS A LIMITED LICENSE UNDER PATENTS OWNED BY THE SALK INSTITUTE FOR BIOLOGICAL STUDIES THAT PERMITS IMANIS LIFE SCIENCES TO SELL PRODUCTS CONTAINING WPRE FOR RESEARCH USE ONLY AND NOT FOR ANY COMMERCIAL USES. EXCLUDED COMMERCIAL USES INCLUDE WITHOUT LIMITATION MANUFACTURING, PROVIDING A SERVICE, THERAPEUTIC, DIAGNOSTIC AND PROPHYLACTIC USES, AND ANY OTHER COMMERCIAL USES. USE OF THIS PRODUCT BY A PURCHASER FOR ANY PURPOSE OTHER THAN FOR RESEARCH IS UNAUTHORIZED AND PROHIBITED.

THE SALK INSTITUTE ACTIVELY LICENSES ITS PATENTS FOR COMMERCIAL USES, AND A COMMERCIAL USE LICENSE MAY BE AVAILABLE FOR SALK'S WPRE PATENTS. IF YOU WISH TO INQUIRE ABOUT SUCH A LICENSE, PLEASE CONTACT:

Office of Technology Development
The Salk Institute for Biological Studies
10010 North Torrey Pines Road
La Jolla, CA 92037
Phone: (858) 453-4100 extension 1278
Fax: (858) 546-8093

THE NIS GENE AND TECHNOLOGY IS COVERED UNDER AN EXCLUSIVE LICENSE TO IMANIS LIFE SCIENCES. RESEARCHERS MAY USE THIS PRODUCT FOR RESEARCH USE ONLY AND NO COMMERCIAL USE IS ALLOWED. NO OTHER USE OR TRANSFER OF THIS PRODUCT OR DERIVATIVES IS AUTHORIZED WITHOUT THE PRIOR EXPRESS WRITTEN CONSENT OF IMANIS. WITH RESPECT TO ANY USES OUTSIDE THIS LABEL LICENSE, INCLUDING ANY DIAGNOSTIC, THERAPEUTIC OR PROPHYLACTIC USES, PLEASE CONTACT IMANIS LIFE SCIENCES FOR SUPPLY AND LICENSING INFORMATION.

Imanis Life Sciences
221 1st Ave SW STE 102
Rochester MN 55902
Tel: (507) 218-2559
Email: info@imanislife.com

THE IMANIS CELL LINES ARE NOT INTENDED FOR USE IN HUMANS. CELL LINES TRANSDUED WITH LENTIVIRAL VECTORS ARE CLASSIFIED AS BIOSAFETY LEVEL 2 REAGENTS AND SHOULD BE USED UNDER THE APPROPRIATE BIOSAFETY LEVEL PER INSTITUTIONAL GUIDELINES.

THE PURCHASER AGREES THAT IMANIS MATERIALS DESIGNATED AS BIO-SAFETY LEVEL 2 CONSTITUTE KNOWN PATHOGENS AND THAT OTHER IMANIS MATERIALS NOT SO DESIGNATED AND ANY PROGENY OR MODIFICATION MAY BE PATHOGENIC UNDER CERTAIN CONDITIONS. PURCHASER ASSUMES ALL RISK AND RESPONSIBILITY IN CONNECTION WITH THE RECEIPT, HANDLING, STORAGE, DISPOSAL, TRANSFER AND USE OF THE IMANIS MATERIALS INCLUDING WITHOUT LIMITATION TAKING ALL APPROPRIATE SAFETY AND HANDLING PRECAUTIONS TO MINIMIZE HEALTH OR ENVIRONMENTAL RISK. PURCHASER AGREES THAT ANY ACTIVITY UNDERTAKEN WITH THE IMANIS MATERIALS AND ANY PROGENY OR MODIFICATION WILL BE CONDUCTED IN COMPLIANCE WITH ALL APPLICABLE GUIDELINES, LAWS AND REGULATIONS.

THE IMANIS MATERIAL, ANY OTHER IMANIS PRODUCTS, AND ANY TECHNICAL INFORMATION AND ASSISTANCE PROVIDED BY IMANIS ARE PROVIDED "AS IS", WITHOUT WARRANTIES OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, MANUFACTURE ACCORDING TO cGMP STANDARDS, TYPICALITY, SAFETY, ACCURACY AND NON-INFRINGEMENT.

IN NO EVENT SHALL IMANIS, ITS PARENTS, SUBSIDIARIES, DIRECTORS, OFFICERS, AGENTS, EMPLOYEES, ASSIGNS, SUCCESSORS AND AFFILIATE (COLLECTIVELY "IMANIS INDEMNIFIED PARTIES") BE LIABLE FOR INDIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES OF ANY KIND IN CONNECTION WITH OR ARISING OUT OF THIS AGREEMENT (WHETHER IN CONTRACT, TORT, NEGLIGENCE, STRICT LIABILITY, STATUTE OR OTHERWISE) EVEN IF IMANIS HAS BEEN ADVISED, KNEW OR SHOULD HAVE KNOWN OF THE POSSIBILITY OF SUCH DAMAGES, INCLUDING, BUT NOT LIMITED TO, LOST PROFITS, COST OF CAPITAL, COST OF SUBSTITUTE PRODUCTS OR CLAIMS OF LICENSEE'S CUSTOMERS FOR SUCH DAMAGE. IN NO EVENT SHALL IMANIS' CUMULATIVE LIABILITY EXCEED THE ACTUAL AMOUNTS PAID BY PURCHASER UNDER THIS AGREEMENT FOR THE TWELVE (12) MONTH PERIOD PRECEDING THE DATE OF THE EVENT GIVING RISE TO THE CLAIM. THE PROVISIONS OF THIS SECTION SHALL SURVIVE THE EXPIRATION OR TERMINATION OF THIS AGREEMENT AND SHALL APPLY EVEN IF THE LIMITED REMEDY SPECIFIED IN THIS AGREEMENT IS FOUND TO HAVE FAILED OF ITS ESSENTIAL PURPOSE.

Quality control by: RLV

Quality Assurance by: SPR

Effective Date: 11-Feb-2016