

# LV-Luc2-P2A-Puro

## Product Description

Product Name: LV-Luc2-P2A-Puro  
 Catalog Number: LV012  
 Lot Number: LV-IM53


Reporter gene: Firefly luciferase (Luc2)  
 Selection gene: Puromycin (Puro)

Quantity: 1 mL  
 Titer:  $2.9 \times 10^7$  TU/mL\*  
 Storage media: Serum free media

Shipping: Dry ice  
 Storage: Store at  $\leq -70^\circ\text{C}$  upon receipt. Freeze-thaw cycles will decrease titer.

Shelf life: One year from date of receipt under proper storage conditions.

This is a ready-to-use lentivirus preparation. The virus encodes the firefly luciferase (Luc2) cDNA under control of the spleen focus-forming virus (SFFV) promoter linked to the puromycin resistance gene (Puro) via a P2A cleavage peptide (see below). The lentiviral vectors are self-inactivating (SIN) vectors in which the viral enhancer and promoter have been deleted. Transcription inactivation of the LTR in the SIN provirus increases biosafety by preventing mobilization by replication competent viruses and enables regulated expression of the genes from the internal promoters without *cis*-acting effects of the LTR<sup>1</sup>.

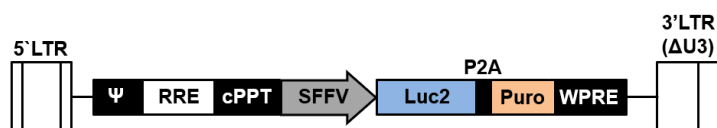

5' LTR: 5' long terminal repeat  
 $\psi$ : RNA packaging signal  
 RRE: Rev response element  
 cPPT: Central polypurine tract  
 SFFV: Spleen focus-forming virus promoter  
 Luc2: Firefly luciferase  
 P2A: P2A cleavage peptide  
 Puro: Puromycin resistance gene  
 WPRE: Woodchuck hepatitis virus posttranscriptional regulatory element  
 3' LTR/ $\Delta$ U3: 3' self-inactivating long terminal repeat

### \*Titration by qPCR:

A WPRE probe-based qPCR assay was used to measure the number of copies of lentiviruses stably integrated into the genome after transduction of HeLa H1 cells (transducing units per mL).

Traditional p24 ELISA titrations measure both functional and non-functional lentivirus particles. However, this method overestimates the functional titer, as the p24 protein pool includes a variable amount of free p24 and p24 associated with non-functional vector particles. This ratio can vary greatly between each lot, so the titration is inherently inaccurate. While qPCR titers may appear lower than p24 ELISA, they are more accurate and functional.

### Safety Precaution:

All culture work with lentiviruses should be performed by trained personnel and performed under BSL2 containment following NIH guidelines.



## Basic Lentivirus Transduction Protocol

Volumes are given for a 6-well plate; increase or decrease as needed. See the Transduction Tips section for additional considerations/modifications.

1. Seed cells in complete medium at an appropriate density to achieve 60-70% confluency the next day (e.g.  $\sim 2.5 \times 10^5$  HeLaH1 cells). Incubate cells overnight in a  $37^\circ\text{C}/5\%$   $\text{CO}_2$  incubator.
2. Thaw lentivirus stock on ice.
3. In a microcentrifuge tube, dilute lentivirus to 1 mL total in serum free media. (See tips below for notes about determining optimal MOIs.)
4. Remove culture medium from cells and replace with prepared lentivirus.
5. Return cells to  $37^\circ\text{C}/5\%$   $\text{CO}_2$  incubator.
6. After 4 hours add 1 mL complete medium to each well and return cells to  $37^\circ\text{C}/5\%$   $\text{CO}_2$  incubator.
7. 3 days after transduction, check transgene expression according to an appropriate protocol. (Note: this lentivirus includes a selection gene; see tips below for details.)

## Transduction Tips

1. To determine the optimal MOI for transductions: plate several wells of the target cells and infect with increasing MOIs (e.g 1, 3, 10, and 30). Typically, primary cells require higher MOIs than established cell lines.
2. Polybrene® (Imanis #REA001) can be added to the transduction mixture to enhance transduction efficiency<sup>2</sup>. The final concentration of Polybrene® in the transduction medium should be 4-8  $\mu\text{g}/\text{mL}$ . Polybrene® can be cytotoxic to some cells and it is not advisable to incubate these cells with Polybrene® overnight; the transduction cocktail may be removed after 3-4 h and replaced with complete media.
3. A spin infection can also be used to increase transduction efficiency<sup>3</sup>. Once the transduction mixture is added to the cells, centrifuge the plate at 800 x g for 30 min at room temperature, before placing the cells in a  $37^\circ\text{C}/5\%$   $\text{CO}_2$  incubator.
4. The presence of serum in the transduction mixture can greatly affect transduction efficiency<sup>4</sup>. In general, lower serum concentrations result in higher transduction efficiencies, though optimizing serum concentrations is recommended for each cell type.
5. The presence of the puromycin resistance gene facilitates selection of transduced cells with puromycin. Selection with puromycin can be performed before or after transgene testing. The appropriate concentration of puromycin to use for selection varies with each cell line and can be determined by performing a kill curve on parental in parallel with transduced cells.

## References

- <sup>1</sup>Miyoshi et al. J Virol. 1998. 72:8150-8157.
- <sup>2</sup>Konopka et al. J Gen Virol 1991. 72: 2685-2696.
- <sup>3</sup>O'Doherty et al. J Virol. 2000. 74:10074-10080.
- <sup>4</sup>Andreadis and Palsson. Human Gene Ther. 1997. 8:285-291.

# LV-Luc2-P2A-Puro




## Certificate of Analysis

Testing performed by Imanis Life Sciences

Test description	Result
Virus titer	$2.9 \times 10^7$ TU/mL
Sterility	No contamination detected
Puromycin selection	Pass QC
Luciferase expression	Pass QC
Endotoxin testing	Pass QC

## Luciferase Expression:

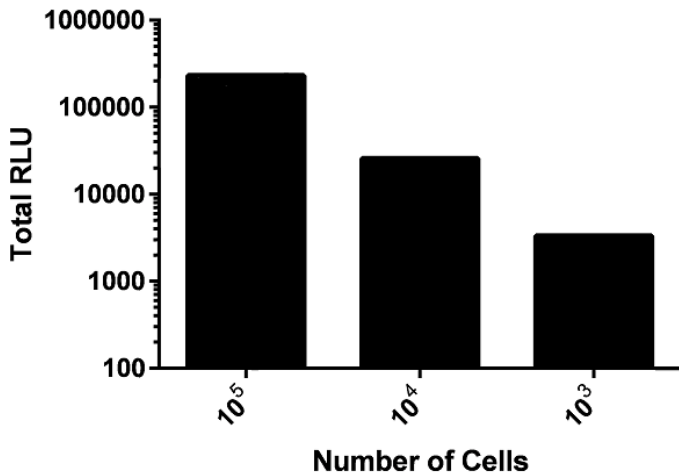

HT1080 cells were transduced with LV-Luc2-P2A-Puro at a MOI of 10. To measure luciferase activity, the transduced cells were seeded into wells of a 96-well plate, 3 mg/mL of d-luciferin was added, and relative light units were measured on a plate reader.

## Legal Disclaimers

### Limited Product Warranty

This warranty limits our liability to replacement of this product. No other warranties of any kind, express or implied, including, without limitation, implied warranties of merchantability or fitness for a particular purpose, are provided by Imanis. Imanis shall have no liability for any direct, indirect, consequential, or incidental damages arising out of the use, the results of use, or the inability to use this product.

For *in vitro* use only. This certificate is a declaration of analysis at the time of manufacture.

### PURCHASER NOTIFICATION

#### LIMITED LICENSE NOTICE – RESEARCH USE ONLY

IMANIS LIFE SCIENCES HAS A LIMITED LICENSE UNDER PATENTS OWNED BY THE SALK INSTITUTE FOR BIOLOGICAL STUDIES THAT PERMITS IMANIS LIFE SCIENCES TO SELL PRODUCTS CONTAINING WPRE FOR RESEARCH USE ONLY AND NOT FOR ANY COMMERCIAL USES. EXCLUDED COMMERCIAL USES INCLUDE WITHOUT LIMITATION MANUFACTURING, PROVIDING A SERVICE, THERAPEUTIC, DIAGNOSTIC AND PROPHYLACTIC USES, AND ANY OTHER COMMERCIAL USES. USE OF THIS PRODUCT BY A PURCHASER FOR ANY PURPOSE OTHER THAN FOR RESEARCH IS UNAUTHORIZED AND PROHIBITED.

The Salk Institute actively licenses its patents for commercial uses, and a commercial use license may be available for Salk's WPRE patents. If you wish to inquire about such a license, please contact:

Office of Technology Development  
The Salk Institute for Biological Studies  
10010 North Torrey Pines Road  
La Jolla, CA 92037  
Phone: (858) 453-4100 extension 1278  
Fax: (858) 546-8093

Quality control by: JDR

Quality Assurance by: RLV

Effective Date: 8-Nov-2019